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



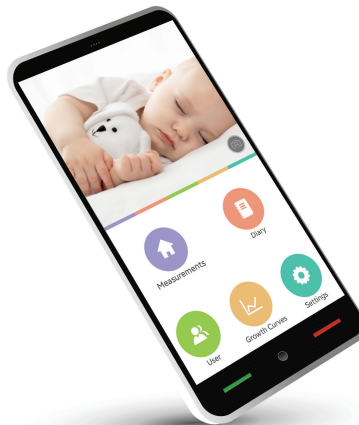
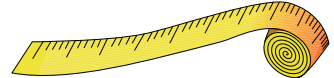
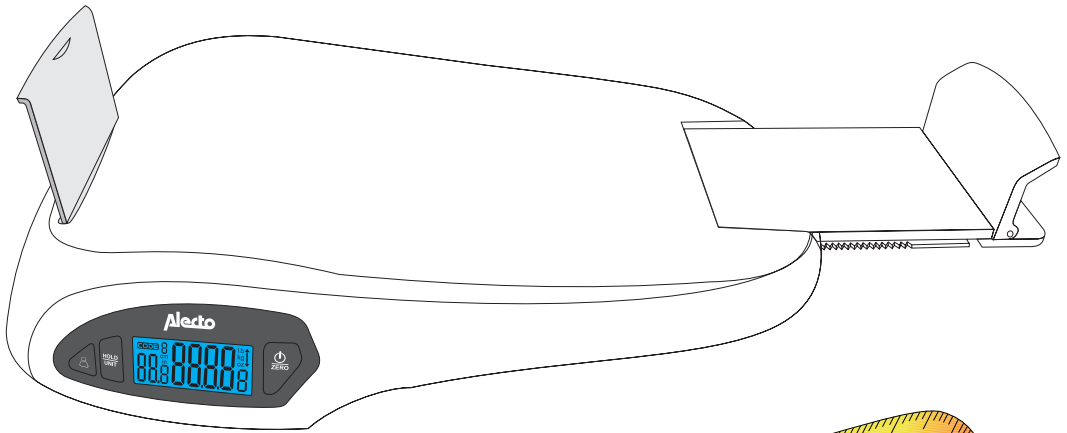
OVERVIEW

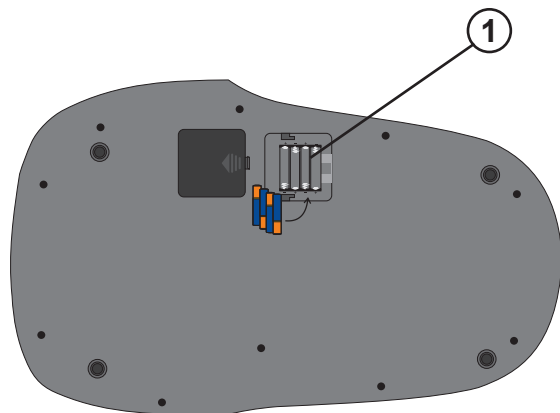1. Battery compartment
2. Main support
3. Feet support
4. Baby selection button
5. Hold/unit button
6. Display
7. Power and 0 reset button
8. Scales
9. Baby memory indicator
10. Baby number
11. Measured weight display
12. Unit of measurement display in Kg
13. Weight progress since previous measurement
14. Measured height display
15. "Ready" icon
16. Measuring tape for measuring the head's circumference

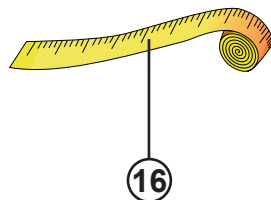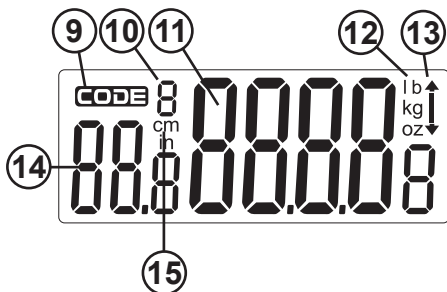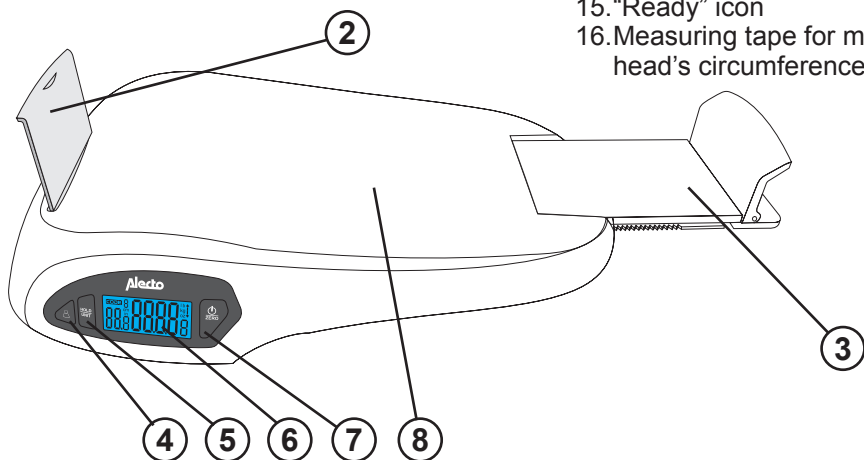INTRODUCTION

The Alecto BC-60BT baby scales allow you to quickly and easily determine the weight and height of your baby. In this way you can't just keep track of how fast your baby is growing, you can also determine how much food your baby has consumed by weighing your baby before and after feeding time. You may use the free app to save all measurements on an Android or iOS device. You may also input the head's circumference measurement. You may then keep track of the growth curve and compare this with the WHO (World Health Organization) growth curve from the birth up to 2 years old.

SAFE USE


- Never leave your baby unattended during measurements.
- The baby scales are not a toy. Never leave your baby alone with the scales; store the scales when they're no longer needed.
- The baby scales are a precision instrument. Never drop the scales and prevent them from becoming wet.

INSTALLATION

Power supply:

The BC-60BT is powered by batteries:

1. Put a towel (to prevent any damage) on a table and then put the scales on the towel upside down.
2. Open the battery compartment and insert 4 batteries, 1.5V, size AAA, Micro or LR03, according to the markings inside the battery compartment. Please keep in mind the correct polarity (+ and -).
3. Close the battery compartment and turn the scales to put them down in the normal way.

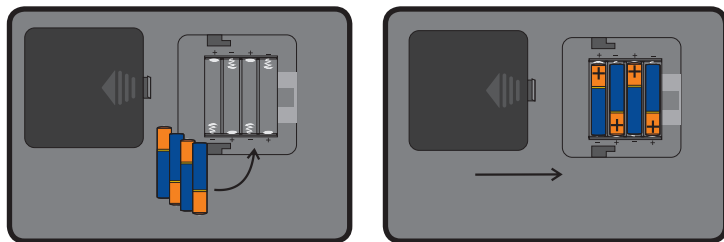The scales are now ready for use.

Always install fully charged batteries of the same brand and type, preferably alkaline batteries. Rechargeable batteries provide a lower voltage and are therefore not applicable.


Positioning:

- Carefully place the scales on a sturdy, stable, flat and dry surface.
- For a reliable measurement, the scales must be put down horizontally.
- Ensure that all 4 feet of the scales make contact with the surface onto which you've placed the scales (the scales may not wobble).
- To prevent falling off, never place the scales too close to the edge of the table or cabinet.

USE


The scales must be prepared to be able to measure your baby's height.

Remove the rubber spacer at the beginning of the weighing platform.

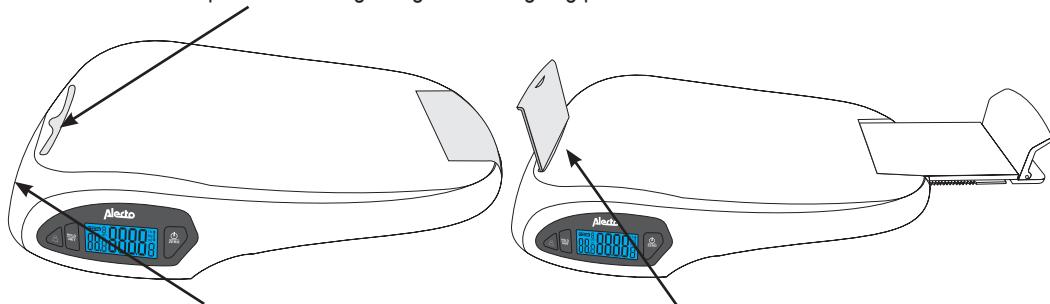Extract the head support from the front of the scales and insert it into the opening you made earlier by removing the rubber spacer.

Flip up the feet support and slide it to the right until there's enough space for the baby to lay on the weighing platform.

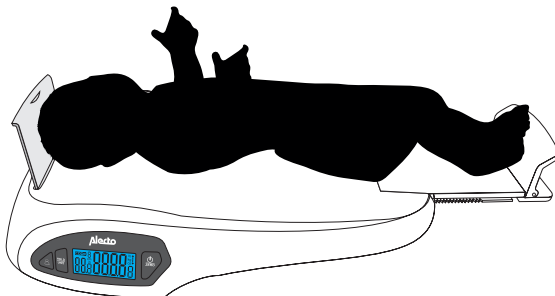Operation:

To prevent undesired enabling or disabling of any functions, you must firmly press down on the buttons. You'll feel a click as soon as the button is pressed down properly.

To switch on:

1. Press and hold the $\frac{\phi}{\text{ZERO}}$ button for at least a second; all digits of the scales' display will now illuminate.
2. The scales are ready for use as soon as the 0.000 kg digits illuminate.

To weigh:


1. Switch on the scales according to the above instructions.
2. Slide the height measure to a height exceeding the height of your baby.
3. Put down the baby in the middle of the scales.
4. As soon as the scales have determined the weight, the weight display will remain fixed.
5. Press the $\frac{\text{HOLD}}{\text{UNIT}}$ button to save the weight; the display flashes twice.

0 reset:

When you first switch on the scales and then cover them with a towel or cloth, the scales will indicate the weight of what you've put down. You can then briefly press the 0 reset button to set the read-out back to 0.000.

To measure height:

1. Switch on the scales according to the instructions above.
2. Slide the height measure to a height exceeding the height of your baby.
3. Put the baby down in the middle of the scales with its head against the head support.
4. Slide the feet support against the baby's feet, making sure to keep the baby lying down as straight as possible.
5. The display will now indicate the baby's height in cm.

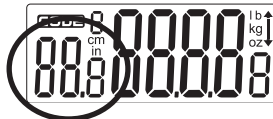To switch off:

If the scales don't determine any change in weight for approx. 60 seconds, the device will switch itself off automatically. You may also press and hold the  button for at least a second to switch off the scales beforehand.

Operation with the app:

Download the BabyGrowing app via App Store (iOS) or the BabyGrowing app via Google Play (Android) and install it on your smartphone or tablet.

The app only works on iOS and Android. Only 1 device may be connected with the scales at any one time. Upon starting the app, you'll be asked to put it into contact with the scales via *Bluetooth®*. You'll need to authorize this.

Baby

Opening the app for the first time.

	Name
	Gender
	Date of birth
	Theme
	User Number

	Joris
	Gender
	Date of birth
	Theme
	

 Baby Girl

 Baby Boy

	Joris
	Baby Boy
	Date of birth
	Theme
	

X	Day	Month	Year	✓
	26	09	2014	
	27	10	2015	
	28	11	2016	
	29	12	1997	
	30	01	1998	

Input the data of the first baby.


Name: Enter the name here.

Gender: Enter the gender here; "Baby Girl" is a girl, "Baby Boy" is a boy. Confirm your choice by pressing the check mark in the gray bar or by pressing Confirm.

Date of birth: Enter the date of birth here. Confirm your choice by pressing the check mark in the gray bar.

Theme: Here you determine the colour of the app when selecting this baby. Choose from purple, green and orange.

User Number: Number of the user. This number corresponds with the **CODE** indication in the display.

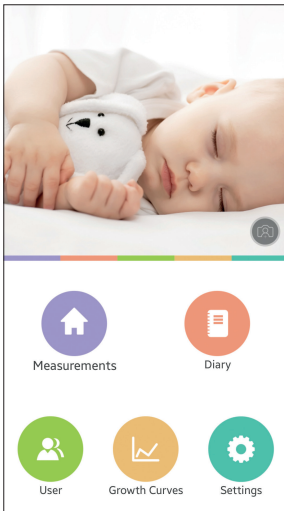
After setting up the first baby, each time you start the app again the alongside image will appear in the start-up screen.

To weigh with the app:

When the data of the baby are entered, press **“Measurements”**, when multiple babies are entered, press **“User”** and select the desired child.

The measurements screen will now appear. Make sure that the **Bluetooth®** function of your smartphone or tablet is activated. Switch on the scales. Put your baby onto the scales according to the instructions in **-To weigh-**. The weight of your baby will now be shown in the weight curve in the measurements screen. Follow the procedure **-To measure height-** to determine the height.

The height is indicated by **“Height”** in the screen. Use the supplied measuring tape to measure the head’s circumference. Slide the triangle to enter the value at **“Head size”**. When all data is known and entered, press **“Save”** to save all data (if the baby is not lying completely still, a notification may appear informing you that the weight isn’t stable and cannot be saved. In this case, press the **“HOLD UNIT”** button on the scales to fix the weight and then press **“Save”** in the app). Press the **“HOLD UNIT”** button on the scales to fix the weight; the display flashes twice. The last value is now saved in the scales. After the display stops flashing, the scales will show the progress since the last measurement. The arrow allows you to quickly determine whether the weight has increased or decreased. When performing multiple daily measurements, the last measurement is saved in the memory.

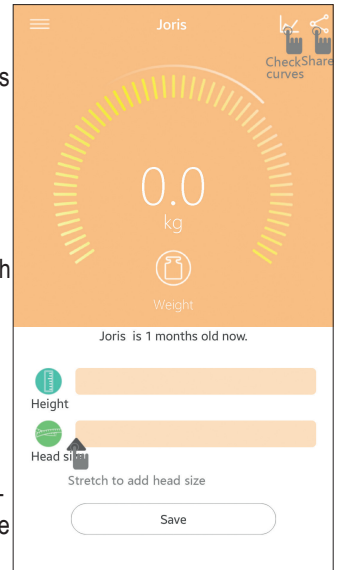To share measurements:

Press the button to share measurements via social media.

To check measurements in the growth curve:

Press the button in the measurements screen or press **“Growth Curves”** in the start-up screen.

You may check the following curves:

Weight/Age - Weight in Kg / Age.

Height/Age - Height in cm / Age.


Head size/Age – Head size in cm / Age.

Weight/Height - Weight in Kg / Height in cm.

Slide the screen upwards to check other curves.

Press the picture of the curve to enlarge the curve.

Press **“All data”** in the upper right corner to adjust, add or remove data.

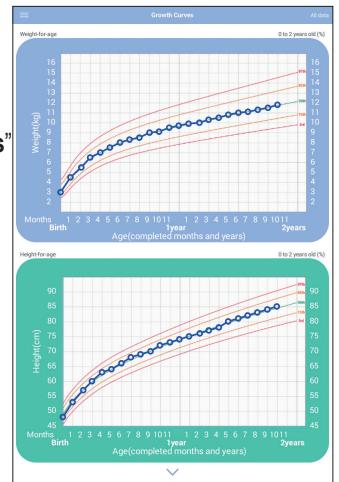To keep a diary:

The app also allows you to keep a diary. Please make sure to select the desired child.

Press the **“Diary”** button to open the diary of the desired child.

Press **“Write”** and **“Entry”** to type a story.

Press the **“Date”** button to choose the desired date. The calendar appears, set and confirm with the check mark.

The story may be entered under **“Note”**.

You may press **+** to make or choose a photo from the file to save together with the story.

To read the diary:

Please make sure to select the desired child.

Press the **“Diary”** button to open the diary of the desired child. Choose the month you'd like to read. Now choose the date you'd like to read or edit. Press **“Done”** after editing to save an adjustment.

Scales memory:

The weight measured last is saved in the memory for each child. Use the baby selection button **“CODE 1...8”** to choose the desired child. When **“CODE”** appears, the scales are set to general mode (no memory). Weight the child and press **“HOLD ONIT”**. The value will now flash several times, the current weight then flashes and finally the previous weight is shown together with the progress icon arrow. The weight has either increased, remained the same or has decreased.

To adjust the start-up photo:

Upon starting the app, the overview menu appears. Above this menu the picture of a baby is shown. You may change this picture in a photo of your baby, child or another photo. Press the  button in the photo. Select whether you want to make a photo or to choose a photo from the file of the smartphone or tablet. The standard baby photo is now changed.

GENERAL USER TIPS AND WARNINGS

App

- The description of the app is merely general and may differ somewhat in the iOS or Android version.

Cleaning:

- Switch off the scales before cleaning them.
- Clean the scales with a moist cloth. Don't use abrasive or chemical cleaning agents. However, you may additionally disinfect the scales with baby-friendly disinfecting agents or soap.
- NEVER immerse the scales into water.

Overload:

- The scales have a maximum weight capacity of 25Kg.
- As soon as these 25 Kg are exceeded, **[Err]** will appear in the display.
- To prevent any damage to the weighing mechanism, never place any objects heavier than 30Kg onto the scales.

Battery power:

- As soon as [Lo] illuminates in the display, the batteries should be replaced.
- Remove the batteries when not using the scales for a long period of time (e.g. several months) to prevent damage caused by battery leakage.

Environment:

- At the end of its lifespan, this product should not be discarded as normal domestic waste, but is must be taken to a collection point for the recycling of electric and electronic equipment.
- Don't throw away exhausted batteries, take them to your local depot for small chemical waste.

Storage:

If the scales must be stored for a long period of time:

- Remove the batteries to prevent damage caused by battery leakage.
- Store the scales horizontally, i.e. NOT on its side.
- DO NOT place heavy objects onto the scales.
- Store the scales in a dry space, within temperatures between 10°C and 60°C.

Service:

- The scales are calibrated upon leaving the factory. The scales don't contain any adjustment controls or spare parts.
- First replace the batteries when you suspect a fault in the device. If this doesn't resolve the problem, please refer to the Alecto site: www.alecto.nl.

Reset:

The scales may be reset to the initial factory settings. All data in the scales will be removed, including the data of the set baby. Make sure the scales are turned off.

- Press and hold the “” button and then press and hold “”. The scales will now turn on.
- Release the “” button and immediately press the “”. Press and hold the “” button.
- The display will show “Clr.” to indicate that the scales are reset.

SPECIFICATIONS

Power:.....4 x 1.5Volt, format AAA, Micro of LR03

Maxim capacity:.....25Kg / division 5 gram

Height measurement range:.....46~80 cm / division 1mm


Memory:.....8 users supported

Support Android and iOS app.

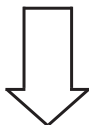The BC-60BT is only suitable for personal purposes

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